


## Separation of Barbiturates

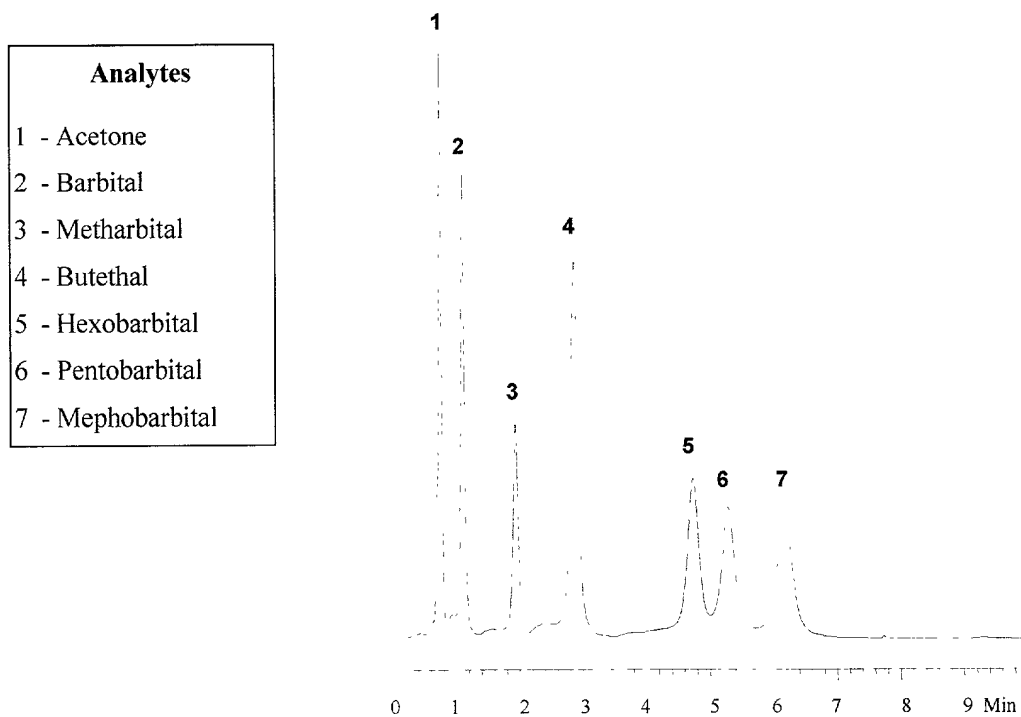

### LC Conditions

Column:	50 x 4.6 mm	Temperature:	30 °C
Mobile Phase: 20/80 A/B		Injection volume:	5 µL
A: ACN		Detection:	254 nm
B: 20 mM Ammonium phosphate, pH 7.0			
Flow rate:	1.0 mL/min.		

Figure 1. Separation of barbiturates using a 50x4.6 mm HPLC column packed with SP-2.

Figure 1

## Separation of PTH-Amino Acids


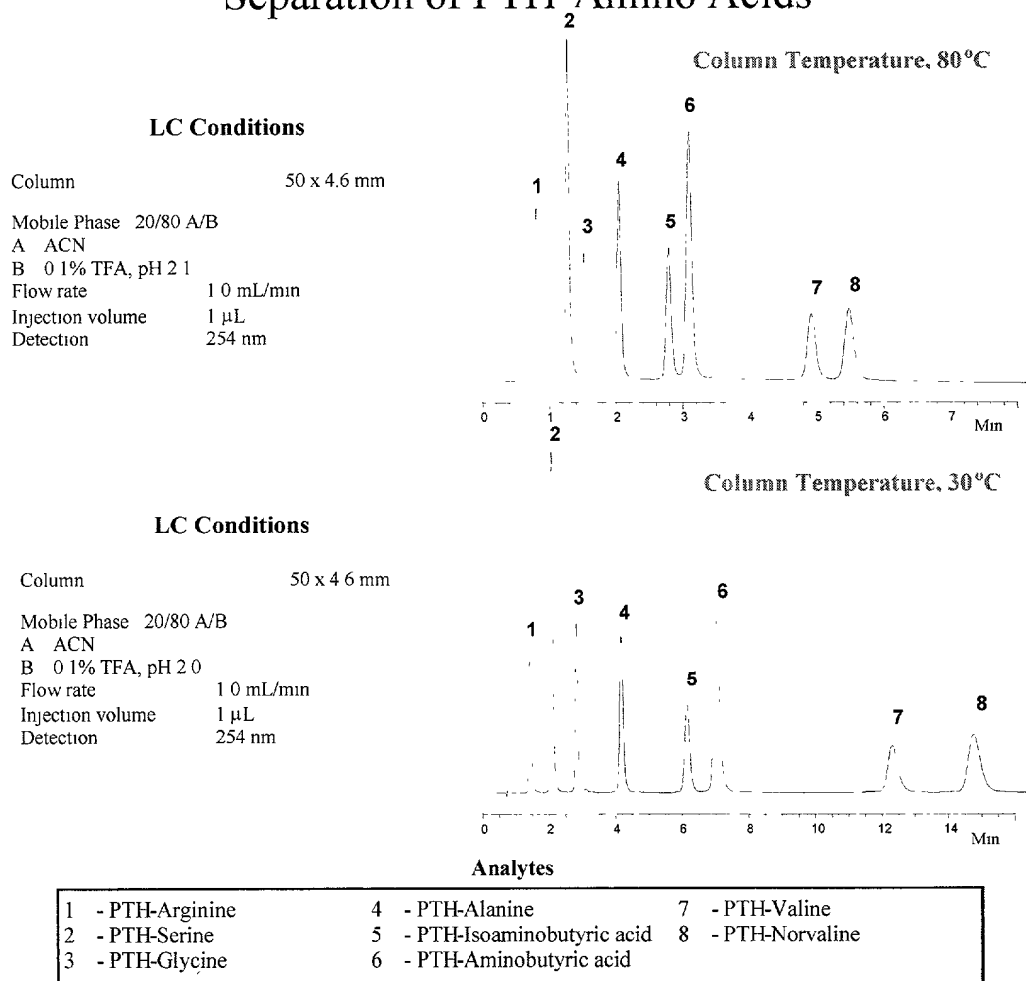


Figure 2. Separation of PTH Aminoacids using a 50x4.6 mm HPLC column packed with SP-2.

Figure 2


## Separation of NSAID's

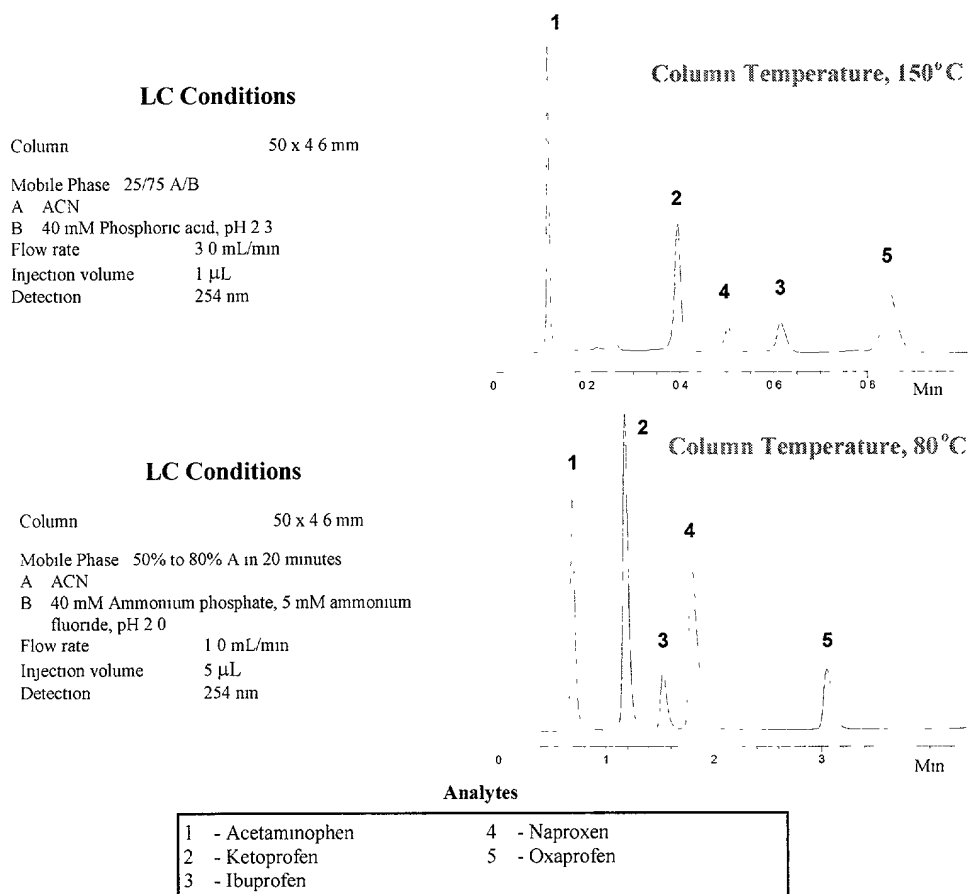

**Figure 3.** Separation of non steroidal anti-inflammatory drugs using a 50x4.6 mm HPLC column packed with SP-2.

Figure 3


## Separation of Ethylbenzene and p-Xylene

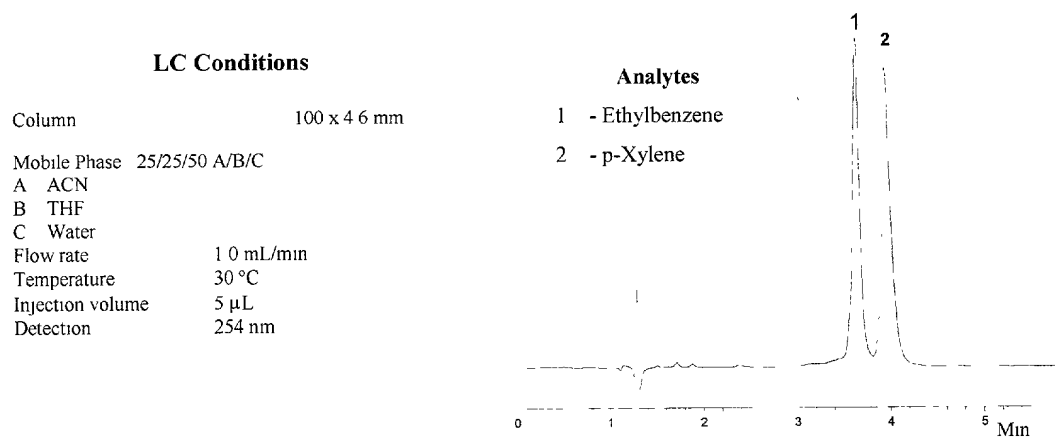

**Figure 4.** Separation of p-xylene and ethylbenzene using a 100x4.6 mm HPLC column packed with SP-2.

Figure 4

## Separation of Beta-Blockers

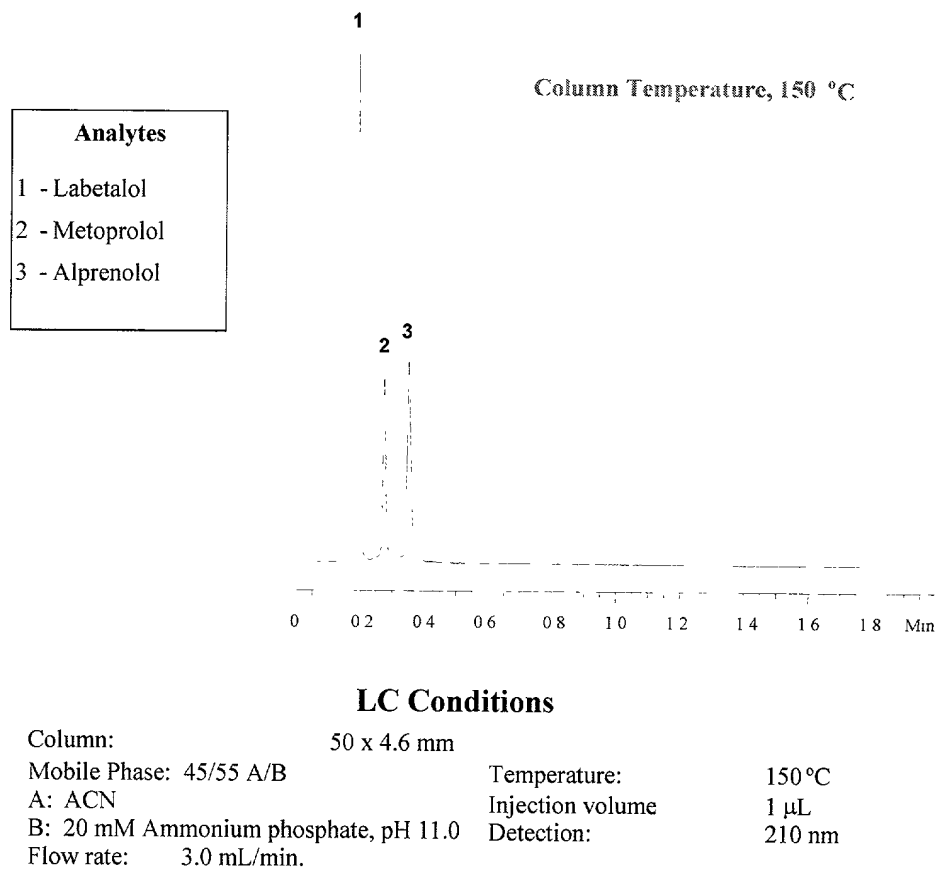

**Figure 5.** Separation of beta-blockers using a 50 x 4.6 mm HPLC column packed with SP-2.

Figure 5